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### **MXC MINI-BOX GENERAL DESCRIPTION**

The Mini-box is made of a galvanised steel plate case painted with highly resistant epoxy paints equipped with a 255 Watt induction motor with internal air recycle through a filter. Internally it is provided with a wedge adjustable in height and width and a PVC tank to collect the residual products.

The upper part is closed by a fireproof laminated plastic plane with rubber edges and in the space below there are 2 slanted and adjustable straticolor arm rests.

The unit is frontally closed by a tempered crystal glass framed by aluminium painted as requested by the client: this crystal is protected by another glass layer.

The mini-box can be located on any working surface. The mobile unit is lighted by two electronically-activated 23 watt screw lamps.

The supplied filter has the following features: 9,200 cm2 nominal filter surface, 20 micron filter nominal capacity.

The panel comprises the following switches:

1 lamp switch 16A

2 suction switch 16A

The mini-box is powered by a 220Volt/50Hz voltage. The electrical input is 46 watt for the lightening and 255 watt for the air suction.

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*The MINI-BOX comes with a one-year period warranty which does not cover any damage or misuse caused by the user.*

*Use the MINI-BOX carefully and in a proper way, avoiding any improper or faulty installation breaching this warranty or safety rules in compliance with its manufacturing standards.*

*Follow this manual very carefully, concerning the proper use and maintenance of the MINI-BOX.*

*Before using the appliance, any operator must be perfectly trained on its functions and technical specifications.*

*The appliance comes with a filter; the package also contains:*


*1 spare filter*

*1 user and maintenance manual and CE certification.*

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*In the mini-box power plug there are two fuses. If the unit does not turn properly on, please check their effective functioning and replace them if needed.*

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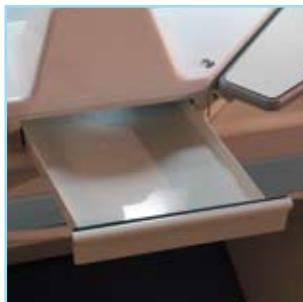*The mini-box is provided with a frontal pull-out drawer. This is useful to collect any surplus waste material. Do not forget to dispose it by simply extracting the drawer from its case.*

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*The internal wedge is adjustable in height operating the mounting screws and in slant pressing it firmly.*

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*To replace or to clean the frontal protection glass, remove the glass holder (see picture on the left). Handle carefully the glass with your right hand, as it may slip when the necessary holder is removed. Finally extract the glass (see picture on the right). To reassemble, repeat the steps in reverse order.*

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*The filter located under the suction unit must be periodically cleaned, depending on the working hours, but at least every 500 working hours. To perform the cleaning, make sure that the power plug is removed from the outlet and the central wedge is lifted. Then undo the hooks on its sides, extract the support (see picture) and remove the filter (see picture) and clean or replace it. The filter must be replaced after 1 year of use. Moreover, make sure that the suction opening is free from any waste material to allow a proper functioning.*

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*The wrist rest slant is adjustable by turning the knob located below it (see picture).*

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**IRIDE INTERNATIONAL S.R.L. Amministratore unico, sig. Enzo Tonin.**

## DECLARATION OF CONFORMITY

We, *IRIDE INTERNATIONAL S.r.l. – furniture for dentist offices and dental laboratories, with registered headquarters in Via San Giuseppe 5 - 33082 Azzano Decimo - Località Tiezzo - PN:*

**HEREBY DECLARE UNDER OUR OWN RESPONSIBILITY THAT THE PRODUCT “MXC MINI-BOX” TO WHICH THIS DECLARATION RELATES WERE DESIGNED AND BUILT IN EXCLUSIVE MANUFACTURING IN CONFORMITY WITH THE FOLLOWING HARMONISED STANDARDS:**

**IEC EN 60601-1(IEC 65-5) of 1991 and its 1994 amendments**

**IEC EN 60204-1(IEC 44/5 and 44/6)**

**AND THE REQUIREMENTS OF THE DIRECTIVE:**

**89/392 of 14/06/1989 and its amendments** relating to the harmonisation of the member states' laws for machinery; **89/336/EEC** electromagnetic compatibility; **89/7323/EEC** high and low voltage.

Machinery excluded by Annex 4 of **Directive 89/392/EEC** in compliance with the general requirement **EN 392 CE 89/392 EEC, CE 89/336 EEC, CE 89/7323 EEC, EN292**

**IRIDE INTERNATIONAL S.R.L. Governing director, Mr. Enzo Tonin.**

**IRIDE INTERNATIONAL s.r.l.**

Via S. Giuseppe 5 | 33082 Azzano Decimo - PN

tel. +39 0434 633242 | fax +39 0434 632546

welcome@irideinternational.com | www.irideinternational.com